




Thank you for purchasing MIYAKIWA
Massage Chair.
This product has a variety of massage
functions, which can be used to relieve
fatigue and improve body function.

MASSAGE CHAIR

instruction manual

MC-5108

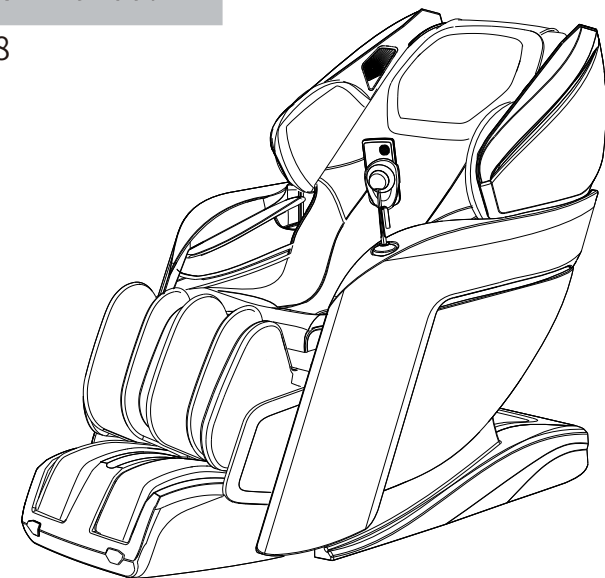

- Please read this manual before use, so as to use the machine correctly.
- Keep this instruction manual in a place convenient for the user to search for, in case of emergency use.
- This machine is only for home use. If it is used for non house hold purposes (such as commercial use), please consult the sales shop or dealer when purchasing this product

RoBo Massage Chair
MC-5108

Manufacturer: MIYAKIWA Health Technology Co., Ltd.

Seller : IIRELAX LIMITED
Unit4, 4 Laurenson Road Hobsonville, Auckland, New Zealand

IRELAX GROUP PTY LTD
Level M/60 Archer St, Chatswood, New South Wales, Australia
2067



Safety Precautions

Please read the safety precautions carefully before use .

Types of illustrations


 Warning	Wrong operation may cause serious injury.
---	---

Types and contents of marks

	Denotes prohibited behavior. (The left figure shows that disassembly is prohibited)
	The content of an instruction that represents a coercive act. (The left figure shows pulling the power plug out of the socket)

After reading, please keep this instruction manual in a place convenient for the user to find, in case of emergency.

Points for attention at the place of use

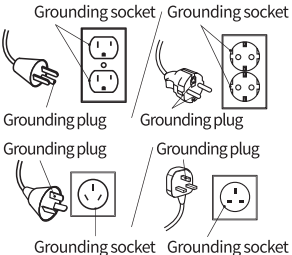 Be careful	Do not use it in the bathroom or sauna and other places with heavy humidity. → Otherwise, electric shock or fault may be caused.
 Force	Place the machine on a flat surface for use. → Otherwise, the massage chair may topple and cause accidents.
 Circuit grounding	- The machine must be grounded. If the machine fails or is damaged, the risk of electric shock can be reduced. A grounding plug is installed in the power cord of the machine, and it can be grounded by inserting the power plug. - The power plug must be plugged into a properly installed and grounded socket completely in accordance with local laws and regulations. If you have questions about whether the machine is properly grounded, you can ask a professional to check. - The grounding plug of this product is shown in the right figure. If it does not match the socket in your home, you need to connect the adapter sold in the market. 

Safety Precautions

Precautions before use

Warning



Be careful

- Do not use a power supply other than AC 220-240V.
→ Otherwise, fire or electric shock may be caused.
- Do not use transformers.
→ Otherwise, failure or electric shock may occur.
- Do not use the machine when the power cord or plug is damaged or the socket is loose.
→ Otherwise, electric shock, short circuit and fire may be caused.



Force

- The following persons must consult a doctor before use.
 - ◎ cardiac ◎ People with osteoporosis, spinal fracture and acute pain
 - ◎ People with severe peripheral vascular circulation disorder caused by diabetes and other factors, resulting in sensory impairment
 - ◎ People with trauma at the treatment site
 - ◎ People whose body temperature is higher than 38 °C (high fever period), for example: the period of severe acute inflammatory symptoms (tiredness, chills, unstable blood pressure, etc.). When the body is weak.
 - ◎ People in the unstable period of early pregnancy or shortly after delivery
 - ◎ People who are using implantable medical devices that are vulnerable to electromagnetic interference such as cardiac pacemakers
 - ◎ People with abnormal spine and bent spine
 - ◎ People who want to use it on the part that has been treated by a doctor or the part of a disease that is being treated
 - ◎ Patients with malignant tumors ◎ People with dizziness, instability and tinnitus
 - ◎ Patients with visceral diseases (gastritis, enteritis, hepatitis) and other acute diseases
 - ◎ People who are being treated by doctors, especially those who feel abnormal
 - ◎ People with rheumatism, gout and other causes of joint deformation
 - ◎ People with inflammatory reaction due to contusion, muscle strain, etc.
 - ◎ People suffering from tenosynovitis or possibly from tenosynovitis
 - ◎ People considered to have lost the sense of temperature
- Before massage, be sure to confirm whether the fabric at the massage part is cracked or damaged. Please check other parts to ensure that the fabric is not damaged. If you find any damage on the fabric, regardless of the size of the damage, please unplug the power plug and contact the dealer who bought this product.
Please never remove the fabric of the massage part or use it when the fabric is damaged, otherwise there will be serious potential safety hazards.
→ If the machine is used when the fabric is damaged, it may cause injury or electric shock.
- The elderly or people with weak bones, even if they do not have specific diseases should consult a doctor before using the massage chair
- Do not let children and pets use or play with the machine.

Safety Precautions

Precautions before use

Be careful



Prohibit

- Do not use the machine in a way not described in this manual. In addition, please do not use this machine with other treatment equipment, or use it after covering with electric blanket
→ Otherwise, it may cause adverse effects on the body or injury.
- Please do not move the massage chair under manned condition.
→ Otherwise, the massage chair may fall over and cause accidents or injuries.
- Do not sit on the chair with the back flat.
→ Otherwise, excessive pressure may be brought to the leg components, resulting in failure. In addition, an accident or injury may occur if the massage chair collapses.
- Do not use anything other than accessories when using this machine.



Force

- Before inserting the power plug into the socket, make sure that the main power switch is in the "O" state.
→ Otherwise, failure or electric shock may occur.
- Please insert properly the power plug and connector
→ Otherwise, electric shock, short circuit and fire may be caused.
- At the beginning, excessive massage stimulation may cause injury. Especially for the elderly or those with weak bones, please pay attention to massage.
→ At the beginning, please start the massage with a light massage program.
- For safety, please check the position of the device. If the device is not in the storage position, please turn on the main power switch, and then slowly sit down after confirming that it has been stored.
- Do not squat or sit on the calf-rest assembly.
→ Massage chair tilting may cause accidents, injuries or failures.
- If the machine has not been used for a long time, please read this instruction manual before using it again to confirm that the machine can operate normally.

Safety Precautions

Precautions in use

Warning



Prohibit

- Do not wind the power cord on the massage chair or step on the body and calf-rest assembly. In addition, please do not damage, damage, process, over bend, stretch or twist the power cord
 - Otherwise, the wire may be damaged, causing fire or electric shock.
- Do not spill water or other liquids on the massage chair.
 - Otherwise, electric shock, short circuit and fire may be caused.
- Do not sit on the massage chair or operate with wet hands when there is water on your body.
 - Otherwise, electric shock or fault may be caused.
- Do not allow children, people unable to express their own will or people with limited mobility to use the massage chair. In addition, do not let children or pets play on or around the massage chair (behind the chair back, under the leg components, etc.), or squat or sit on the massage chair.
 - Otherwise, accidents or injuries may be caused.
- Do not use the massage chair when more than 2 people are sitting on the chair. Especially when holding a child (baby, etc.).
 - Otherwise, accidents or injuries may be caused.
- Stop using if the fabric is damaged, and pull out the power plug from the socket
 - Otherwise, accidents, injuries or electric shock may result. For maintenance, please contact the dealer when purchasing this product
- Do not touch the LED directly when the light is on or just off.
 - Otherwise, electric shock or fault may be caused.



Force

- Before starting the massage and when standing up and lowering the back of the chair, make sure there are no people (especially children) and pets around. In addition, please be careful not to clamp your body or objects behind the chair or under the leg components.
 - Otherwise, accidents or injuries may be caused.
- When your body or object is caught, please press the seat back adjustment button to return the seat back to a safe position.
- Do not cover or touch the LED part with cloth, paper, plastic stool and other flammable articles.
- Do not insert metal objects (hairpin , pins) into the gap of the machine (especially the LED part and power cord).
 - Otherwise, electric shock, short circuit or fire may be caused.

Safety Precautions

Precautions in use

Be careful



Prohibit

- Do not use the massage chair directly with bare skin.
- Do not use the massage chair for purposes other than massage and chair.
 - Otherwise, accidents or injuries may be caused.
- Do not put your hands in the movable parts of the massage wheel, such as the up and down motion parts, the chair back, etc. In addition, do not wear hard objects such as belts when using.
 - Otherwise, it may cause injury or failure.
- Do not use on the head, chest, abdomen, joints (elbows, knees, etc.).
 - Otherwise, it may cause physical discomfort or injury.
- Do not operate the massage chair when nobody is around.
 - Otherwise, accidents or injuries may be caused.
- The calf-rest assembly can only be used for legs.
 - Otherwise, it may cause injury.
- When massaging the leg, do not put the leg in again forcibly after the leg is pulled out of the calf-rest assembly.
 - If forced, it may cause failure.
- The armrest assembly can only be used on the wrist.
 - Otherwise, it may cause injury.
- When massaging the wrist, do not force it into the wrist after the wrist is extruded from the armrest assembly.
 - Forced insertion may cause failure.
- Do not look directly at the LED of the side panel assembly
 - Otherwise it may cause injury.



Force

- In case of emergency or abnormal operation of the massage chair, please press the [Safety Stop] button on the remote control to stop the massage.
- If you feel severe pain or physical abnormalities during use, please stop using immediately and consult a doctor.
- If you feel that the massage is too stimulating, please reduce the irritation according to this manual. If the stimulation is still too strong after adjustment, please stop using and contact the dealer who bought the product (recorded on the back cover).
 - If you continue to use the product even though it is too stimulating, you may get injured.
- Please avoid falling asleep during massage
 - Otherwise it may cause injury.
- The massage time shall not exceed 40 minutes every day. Do not massage the same place for more than 5 minutes.
 - Continuous use for a long time will cause excessive stimulation, resulting in adverse effects or injuries. Please control the total massage time within 60 minutes every day.
- Before use, make sure there are no people (especially children) and pets around.
 - Otherwise, it may cause injury or failure.
- Do not remove the lower back cushion during massage.

Safety Precautions

Precautions in use

⚠ Be careful

 Unplug the power plug	● In case of power failure or possible power failure, please stop using immediately, turn off the power switch and unplug the power plug from the socket. → Otherwise, accidents or injuries may occur when the power supply is restored. ● When thunder may occur, please stop using immediately, turn off the power switch, and unplug the power plug from the socket. ● When the massage chair does not work or feels abnormal, to prevent accidents, please stop using it immediately, turn off the power switch, unplug the power plug from the socket, and contact the dealer who purchased the product for maintenance. → Otherwise, electric shock, electric leakage or short circuit may be caused, which may cause fire.
--	--

Precautions after use and when not in use

⚠ Be careful

 Prohibit	● Do not use the power cord for equipment other than this machine. ● After using this machine, make sure the seat back is completely upright. ● Do not put articles on the massage chair ● Do not replace the LED lamp.
 Force	● After use, be sure to turn off the power switch and unplug the power plug from the socket. → Otherwise, accidents or injuries may be caused by children or pets. ● After use, please clean and place in a dry place.
 Unplug the power plug	● When unplugging the power plug from the socket, do not pull the power cord directly, but hold the power plug and pull it out. → Otherwise, electric shock or short circuit may be caused, which may cause fire. ● Be sure to unplug the power plug from the socket before cleaning and maintenance. Also, do not plug or unplug the power plug with wet hands. ● When not in use, be sure to unplug the power plug from the socket. → Dust or moisture may cause aging of wire insulation, which may cause electric shock, electric leakage or fire.

Safety Precautions

Other points for attention

⚠ Warning

 Prohibit	● Do not stand on the massage chair. Do not give strong impact to the massage chair from above, such as sitting outside the chair seat. → Massage chair tilting may cause accidents, injuries or failures.
 Do not disassemble	● If the massage chair breaks down, please contact the sales shop or dealer and do not repair it by yourself. ● If the power cord is damaged, please entrust the dealer who bought the product to replace it to avoid danger. ● Do not modify the massage chair. → Otherwise, fire, electric shock or injury may be caused.

⚠ Be careful



 Force	● If you transfer or lend the massage chair to others, please be sure to hand over this manual together with the massage chair to the other party.
--	--

Be careful

- When discarding this product and its accessories, please consult the nearest municipal window for the correct method of discarding.

Product Structure

The illustration is a schematic diagram. Different models have different shapes.
Function distribution diagram

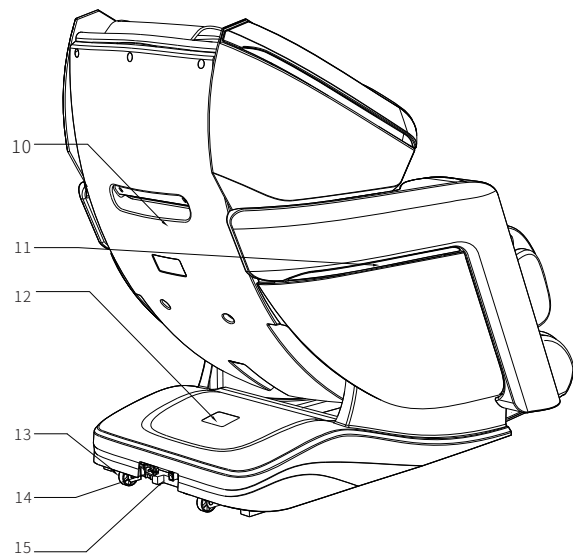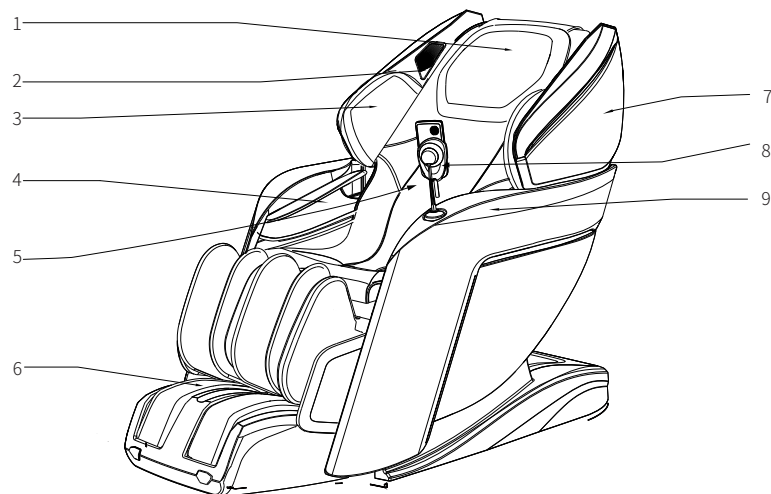

- 1.Pillow
- 2.3D digital audio
- 3.Upper arm airbags
- 4.Arm airbags
- 5.Backrest cushion
- 6.Calfrest
- 7.Space capsule
- 8.5 inch touch
- 9.Armrest
- 10.Backrest cover
- 11.LED light
- 12.Seat frame back cover
- 13.Caster
- 14.Power jack
- 15.Power switch

Operating Instructions


Function

- This product is designed with a set of intelligent massage manipulator which can walk up and down, knead and tap, and silent design.
- It is designed with shoulder position automatic detection and adjustment, it can automatically calculate acupuncture points according to different shoulder heights so as to offer humanistic and scientific massage.
- There are 4 scenario modes, each containing 5 automatic programs, for a total of 20 automatic programs.
- Full Body Mode includes 5 programs: Master's choice, Fatigue recover, Deep Shiatsu, Lymphatic drainage, Elderly's vitality.
- Stretch Mode includes 5 programs: Early morning centripetal, Artisan stretch, Full body stretch, Warm-up exercise, Cooling exercises.
- Soothing Sleep includes 5 programs: Japanese Style Massage, Easy Swing, Sleep aid, Air pressure activation, Swing mode.
- Concentrated Areas includes 5 procedures: Brain relax, Neck&Shoulder, Lower back, Hip, Limb massage.
- The manual interface has four sub-functions: Massage techniques, Mechanism range, Additional function, Airbag massage.
- The lower backrest frame uses electric push rods combined with slide rail technology to maintain zero space and save space to the maximum extent.
- The backrest frame guide rail adopts SL curve guide rail technology, which completely fits the curve of the human back, allowing you to enjoy a close-fitting massage from the head to the thighs.
- Equipped with intelligent voice control function: the left speaker has a built-in voice control microphone, so that you only need to move your mouth, and the massage chair can press your wishes.
- Wireless Bluetooth music function: Built-in 3D digital speakers on the left and right upper arms, you can listen to music and relieve stress while massaging.
- Foot roller massage: effectively relieves body fatigue and promotes sleep.
- Patented calf rubbing: imitating Thai techniques, reaching deep sore muscles and fully relaxing the calf muscles.
- Automatic expansion and contraction of the feet: tailored matching according to different heights to customize an exclusive massage experience.
- Equipped with calf frame and backrest frame lifting functions.
- Equipped with a 5-inch touch screen controller.

Operating Instructions

I. Power on the whole machine

- 1 Schematic diagram of power on
- 2 Schematic diagram of connecting line of the whole machine
- 3 Schematic diagram of power switch of the whole machine

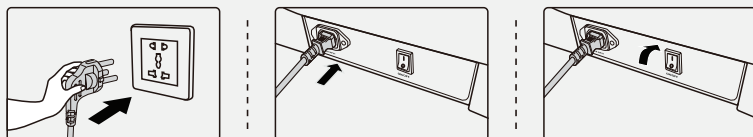

II. Massage function control operation

1. Turn on the machine

- Massage chair power on in standby mode, press the power button
- Display the power-on interface:



Operating Instructions

Home page feature introduction





- | | |
|--|--|
| 1. Pause/Run button: Pause/Continue massage function | 12. Back lift button: Lift the back. |
| 2. Leg lift button: Lift the calf | 13. Air pressure display. |
| 3. Leg lowering button: Lower the calf | 14. Movement position display. |
| 4. Heating part display | 15. Back down button: Lower the back. |
| 5. Strength key of the movement: controls the strength of the movement | 16. Pressure intensity key: controls the pressure intensity. |
| 6. Time adjustment and countdown area | 17. Heating key: controls the heating gear. |
| 7. Hertz movement key: turn on/off movement vibration function | 18. Roller key: controls the force of the foot roller. |
| 8. Menu Bar: Display/Switch Function | 19. Rubbing key: controls the rubbing force of the calf. |
| 9. Settings key: Enter the settings center | 20. Pause/Run button: Pause/Continue massage function. |
| 10. LOGO | 21. Power button: Turn on or off the massage chair. |
| 11. Status bar: Display current mode | |

Note: Hertz movement can only be turned on in automatic mode.

Operating Instructions

2. Automatic lay-down adjustment



- After power on, the touch screen will enter the massage automatic program selection interface, after selecting the automatic program, the massage chair will enter the reclining state and reach the safe angle. then the shoulder position will be detected automatically.

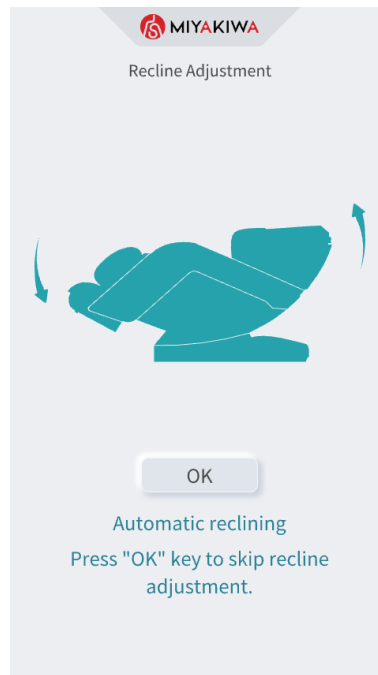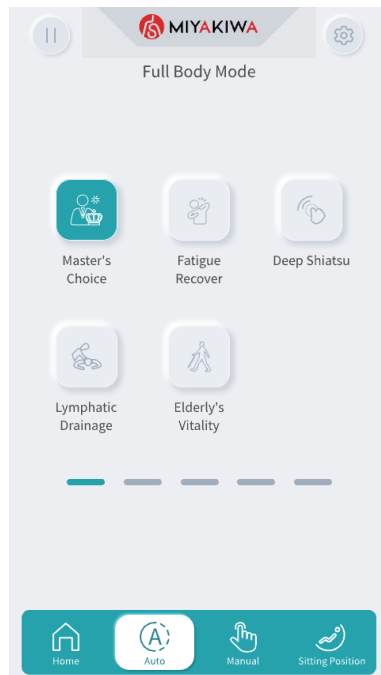

Operating Instructions

3. Shoulder position detection and adjustment and massage position adjustment



- After the manipulator detects the shoulder position, the buzzer prompts "Didi". If it is different from the actual shoulder position, you can adjust it to the appropriate position by clicking "▲" or "▼".
- The shoulder position is adjustable in 11 levels. If there is no shoulder position adjustment within 10 seconds, it will automatically enter the massage function.
- If you press the OK "button", you will directly enter the massage. If the shoulder position cannot be detected, massage will be performed according to the default shoulder position height.

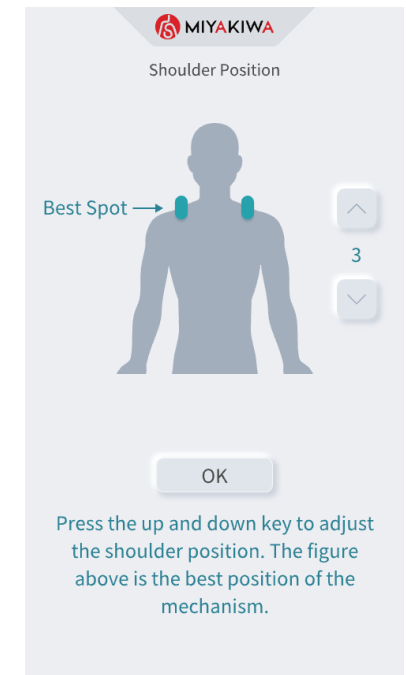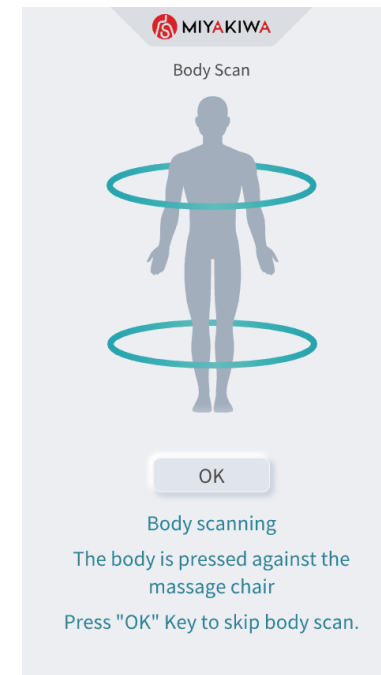

Operating Instructions



4. Automatic program

- Click the "Auto" button to enter the auto program interface, select the auto program icon; there are four scenarios, each scenario contains four auto programs.
- Full Body Mode includes four programs: Master's choice, Fatigue recover, Deep Shiatsu, Lymphatic drainage, Elderly's vitality.
- Stretch Mode includes four programs: Early morning centripetal, Artisan stretch, Full body stretch, Warm-up exercise, Cooling exercises.

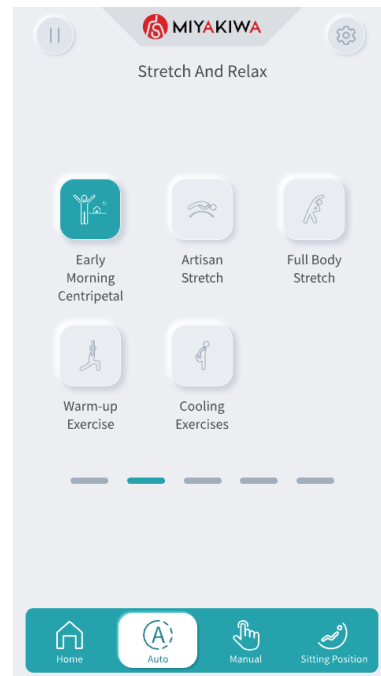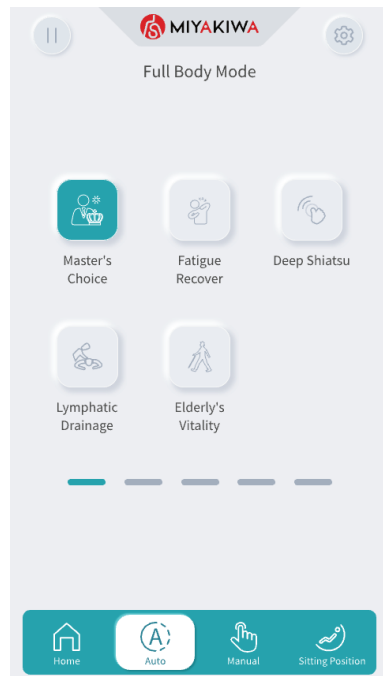

Operating Instructions

- Soothing Sleep includes four programs: Japanese Style Massage, Easy Swing, Sleep aid, Air pressure activation, Swing mode.
- Concentrated Areas includes four procedures: Brain relax, Neck&Shoulder, Lower back , Hip , Limb massage.

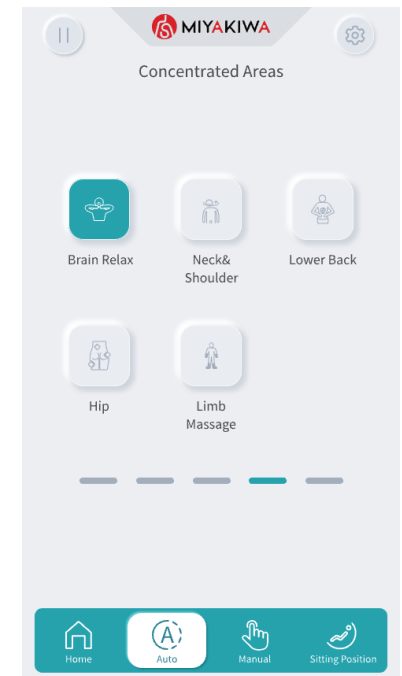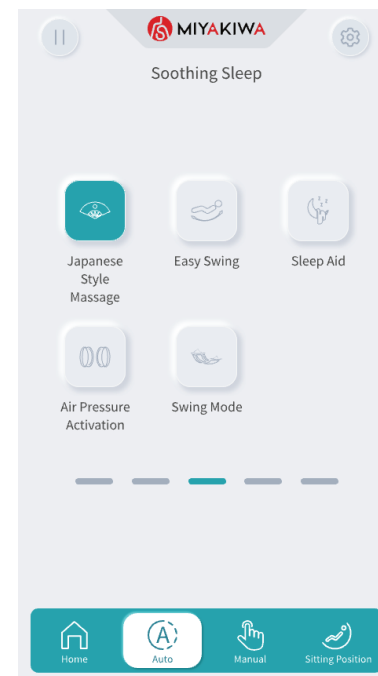

Operating Instructions

Full Body Mode

Master's Choice: A curated combination of massage techniques provides a full-body relaxation experience akin to that from a professional masseur/masseuse.

Fatigue Recover: Targeting muscle soreness after intense exercise or prolonged work, it employs kneading, tapping and percussion techniques to alleviate muscle fatigue.

Deep Shiatsu: Simulating traditional shiatsu therapy, it deeply presses major meridians and acupoints of the body.

Lymphatic Drainage: Using gentle rubbing techniques coordinated with airbag inflation, it promotes lymph circulation and helps the body eliminate metabolic waste.

Elderly's Vitality: Full-body massage promotes blood circulation, bringing a sense of relaxation, comfort and vitality.

Stretch Mode

Early Morning Centripetal: Simulating gentle morning stretches to awaken the body, it mildly extends joints and muscles, activating the body and bringing fresh energy.

Artisan Stretch: A signature stretching program achieves a full-body stretch of legs, hips, waist, and shoulders.

Full Body Stretch: A comprehensive stretching program covers thorough extension from head and neck down to legs and feet.

Warm-up Exercise: For use before exercise, it combines kneading, tapping, and stretching to gently raise muscle temperature and heart rate, as well as to activate joints.

Cooling Exercise: For use after exercise, it helps tense muscles gradually relax, slows heart rate, and reduces delayed onset muscle soreness.

Operating Instructions

Soothing Sleep Mode

Japanese-Style Massage: Employing traditional Japanese massage techniques like shiatsu, rubbing and kneading, it creates a serene and peaceful state of mind, inducing deep relaxation.

Easy Swing: With Zero-Gravity swing mode on, it simulates the comforting sensation of a cradle, calming nerves and releasing stress.

Sleep Aid: With Zero-Gravity mode integrated with extremely gentle kneading, pressing, and stretching, it helps relax the body and mind, guiding towards a sleep state.

Air Pressure Activation: Relying on the airbag system, it performs inflation/deflation compression on various body parts, gently promoting blood circulation in the extremities.

Swing Mode: Under Zero-Gravity swing mode, combined with very light massage and airbag compression, it quickly soothes emotions, promoting deep relaxation and sleep.

Concentrated Areas Mode

Brain Relax: Deep neck massage and vibration relieve head tension, soothe nerves, and alleviate heaviness in the head.

Neck & Shoulder: Focused massage on the shoulder and neck area tackles stiffness, soreness, and knot problems in this region.


Lower Back: Lower back-targeted massage substantially alleviates muscle strain, stiffness, and soreness around the lumbar.

Hip: Aimed at people with sedentary lifestyle, it massages down to the sacrum, tailbone, and buttocks, relieving tension around the pelvis and improving discomfort in the sciatic nerve area.

Limb Massage: "Precisely focusing on the extremities of the arms and legs, it presses via airbags, combined with foot rollers, to vigorously promote blood circulation in the hands and feet."

Operating Instructions



Memory Mode enables recall of 6 preset massage programs. Long press (3s) for storage, single click for activation.

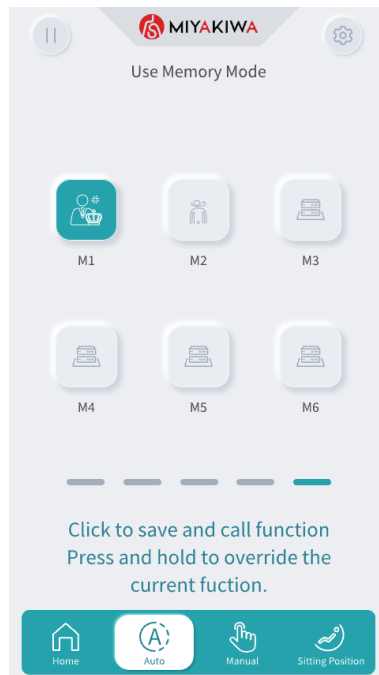

Operating Instructions

5.Manual



- Click the "Manual" button to enter the manual interface. The system defaults to the "Massage Techniques" interface, which can be switched through the slider above.
- There are a total of 9 massage techniques: tapping, kneading, acupressure, percussion, and kneading. There are 5 techniques available for switching in kneading.
- Massage intensity offers 6 adjustable levels. Selection can be made through progress bar clicks or slider manipulation.
- Massage speed provides 5 selectable levels (manual mode exclusive). Slider movement to extreme left terminates function.
- There are three options for the width of the movement. In manual mode, the massage technique can only be adjusted by tapping, tapping, and finger pressure.

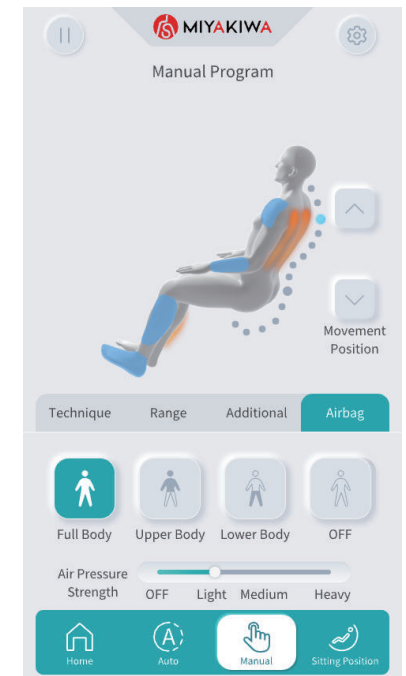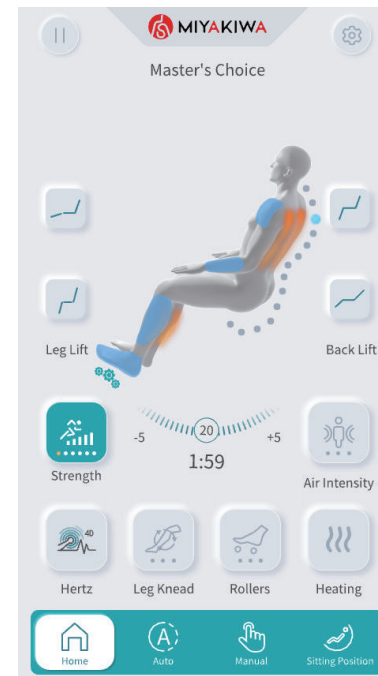

Operating Instructions

6.Range



- Click the "Range" button on the upper side to enter the range interface of the movement, including full back, partial, and fixed point.
- Width: Slide the width adjustment slider to set the massage mechanism to narrow, medium, or wide.

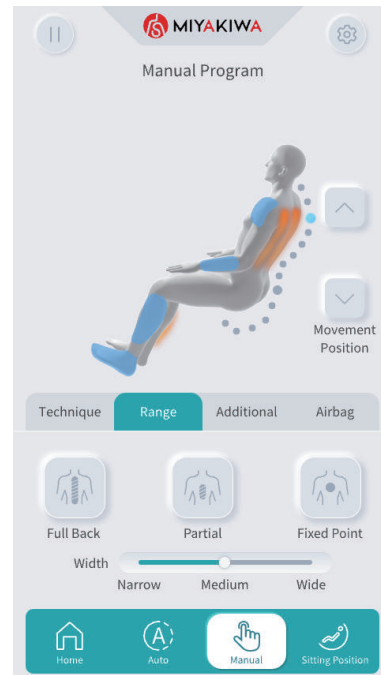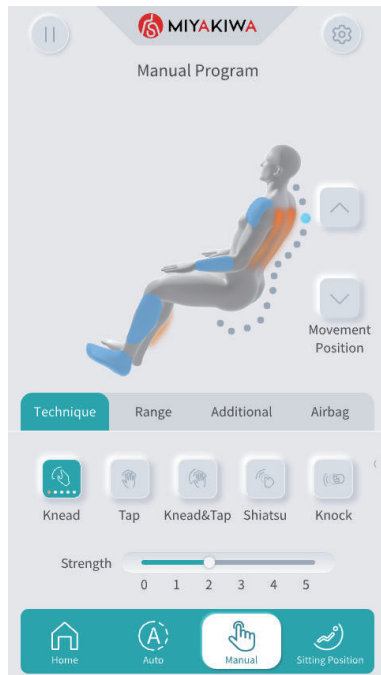

Operating Instructions

7.Additional


- Click the "Additional" button on the upper side to enter the Additional friction interface, where you can adjust the kneading intensity and foot roller speed by clicking the progress bar or dragging the slider.
- There are 3 levels of kneading intensity and plantaroller speed, and the slider is turned off by moving it to the leftmost position.

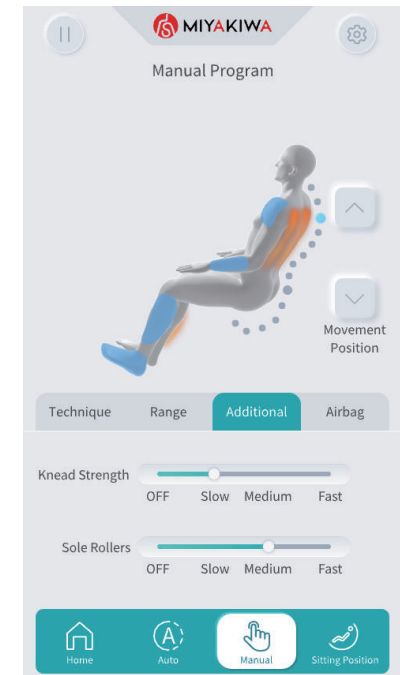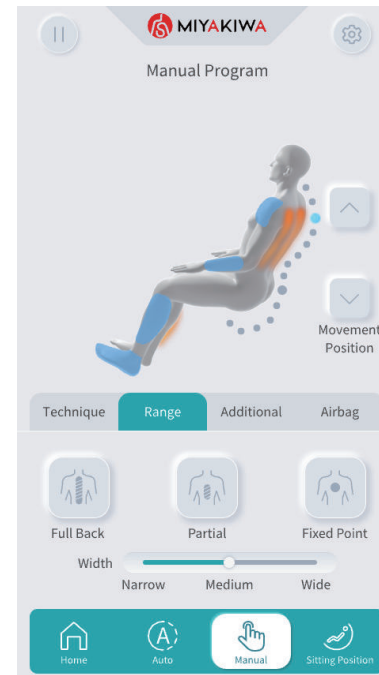

Operating Instructions

8.Airbag

- Click the "Airbag" key to enter the airbag massage interface.
- Airbag Massage: Features whole body, upper body, and lower body modes.
- Air Pressure Intensity: The slider adjusts through three intensity levels. Moving it to the leftmost position turns off air pressure.

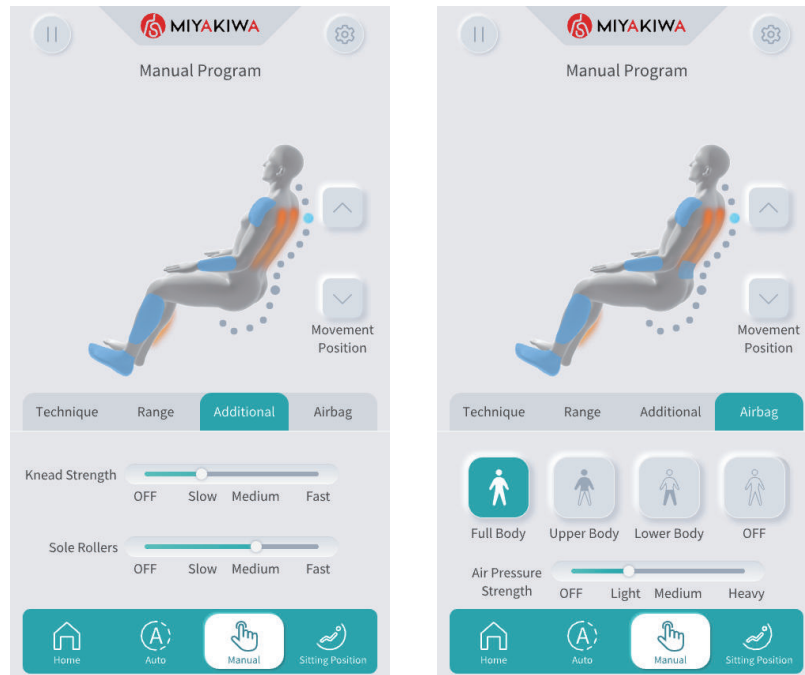


Operating Instructions

9.Seat position adjustment

- Click "Sitting Position " to enter the seat position adjustment interface.
- Long press " " or " "to adjust the leg up .Release the key to stop the leg up.
- Long press " " or " "to adjust the back up .Release the key to stop the back up.
- Long press " " or " "to adjust the leg flexion.Release the key to stop the leg flexion.
- Click "ZEROI "key to enter ZEROI mood or click "ZEROII "key to enter ZEROII mood.Click "ZEROIII " to go to ZEROIII.


Note: When adjusting the lift of the lower leg, the length of the lower leg will be automatically detected after the key release. When the foot retracts against the heel, the foot stops expanding and contracting.(In order to avoid shin frame colliding with the ground when the foot expands and contracts,the shin frame angle must reach a safe angle before the function of expansion and contraction is effective.)

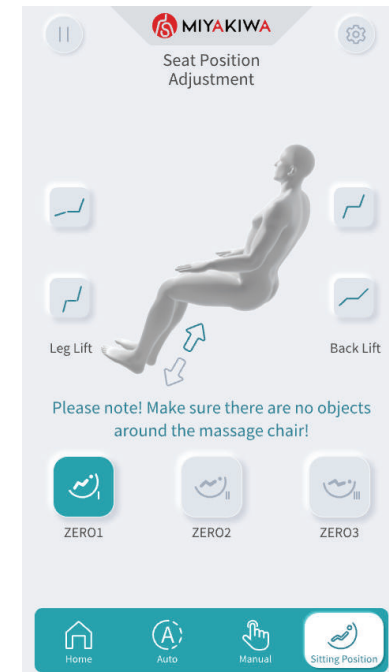

Operating Instructions





10. Automatic program

- Click the "⚙️" button to enter the settings interface.
- Bluetooth: Click the button to enable/disable Bluetooth. After enabling, turn on your phone's Bluetooth, search for the device name "MC-5108", and click to pair. Once connected, you can play music from your phone through the massage chair.
- Shoulder Position Calibration: Click the button to enter the shoulder position calibration interface.
- Elderly Mode: Click to activate Elderly Mode, which simplifies and enlarges the interface, removes complex manual settings while retaining other functions. Click "🏠" in the top-right corner to return to the normal interface

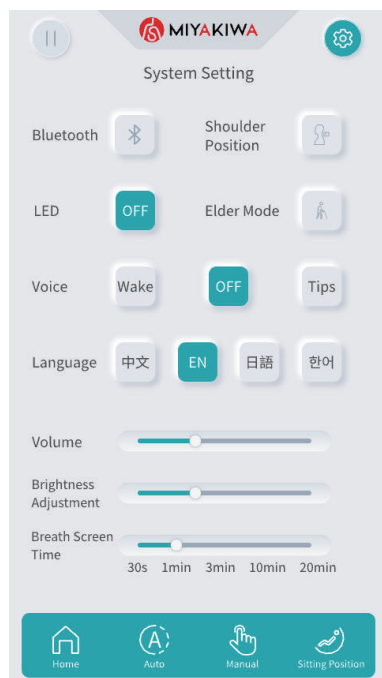



Operating Instructions

Elderly Mode Interface Overview

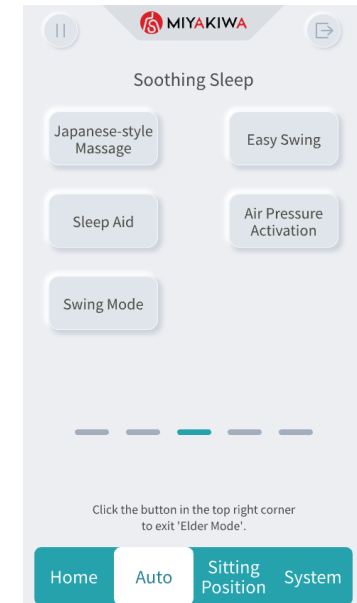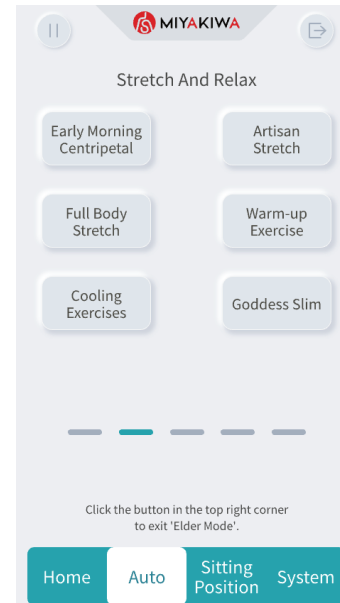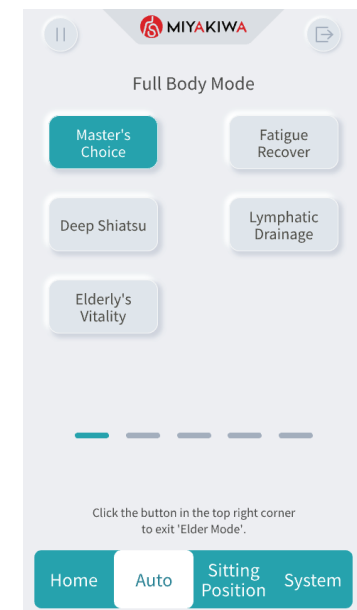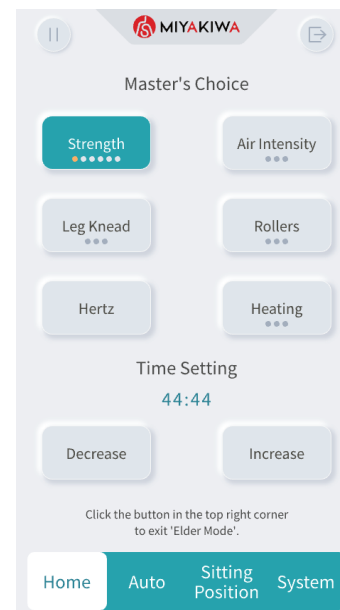


Operating Instructions

Elderly Mode Interface Overview

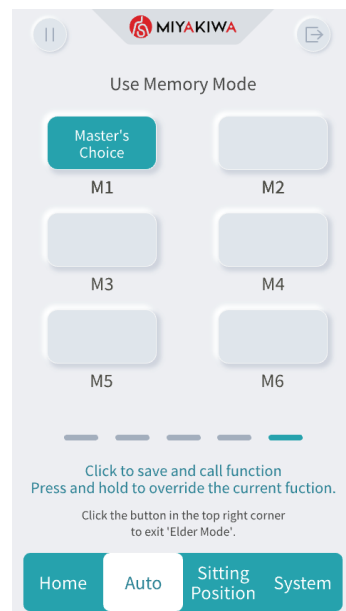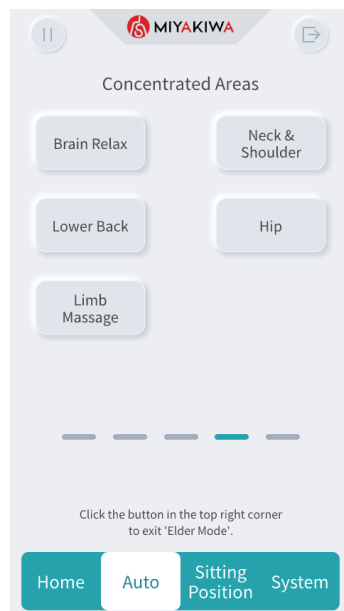


Operating Instructions

Elderly Mode Interface Overview

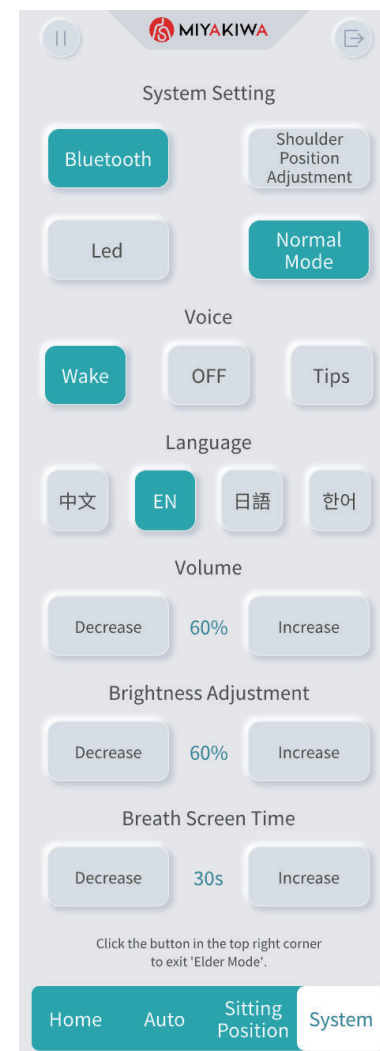

Operating Instructions

- Language: Click to select language: Chinese, English, Japanese, Korean.
- LED: Click the button to turn the ambient light on/off.
- Volume: Click on the progress bar or drag the slider to adjust button volume.
- Brightness: Click on the progress bar or drag the slider to adjust screen brightness.
- Breath screen time: Click on the progress bar or drag the slider to select screen timeout: 30 seconds, 1 minute, 3 minutes, 10 minutes, 20 minutes.

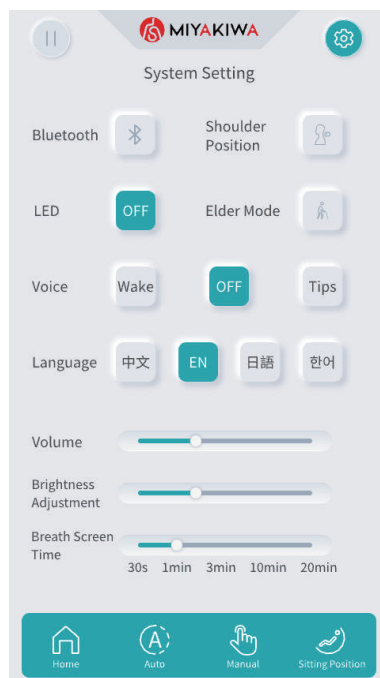


Operating Instructions

- Voice Control Function: After enabling the voice control function, say "Hi what's up" when the massage chair is powered on or in standby mode to activate voice control. Upon successful activation, the chair will respond with "I'm here" and await voice commands within 6 seconds to activate the massage chair.

Note: If no valid voice command is detected within 6 seconds, the voice control will automatically deactivate.

Note: Avoid operating the voice control in noisy environments. When activating the function, use clear Mandarin Chinese at natural speed and strictly follow the predefined voice command phrases.


Note: Voice control enters sleep mode 30 seconds after power-off. Briefly press the power button to wake the function.

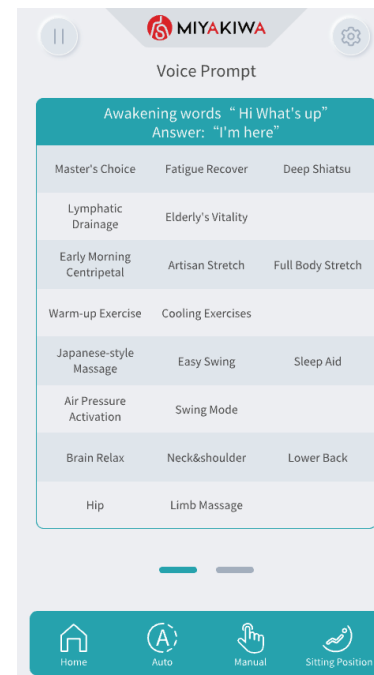

Operating Instructions

Product Installation


1.Unpacking

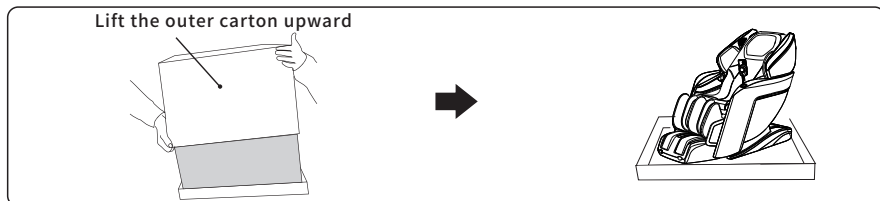

- Cut the packaging straps of the main carton. Lift the outer carton upward to remove it.
- Remove and organize all chair accessories.

2.Moving Method

- Raise the backrest to its highest position.
- Tilt the chair forward by lifting the front end and pressing down on the backrest. Use the rear casters to maneuver. Slowly lower it to the floor at the desired location.

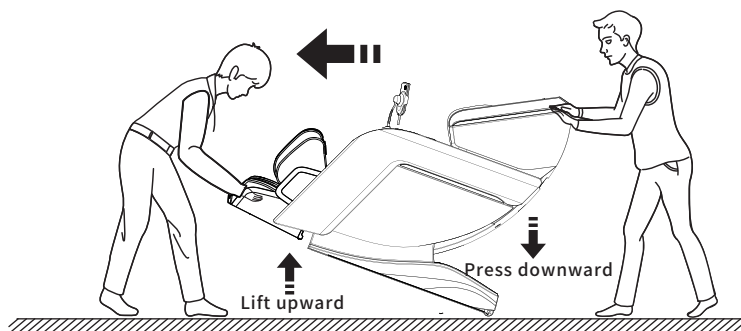

Important:

- Always turn off the power switch before moving.
- Disconnect both the power cord and hand control cable prior to relocation.

Before Use

- This product is packaged with the seat back upside down and automatically resets after power on.
- When the room temperature is low, if you want to increase the indoor temperature and use the massage chair, please gradually increase the temperature. Do not increase the temperature too quickly. When moving the machine from a low-temperature storage place to a warm room, please wait about 1 hour before using it.

→ Otherwise, condensation may occur at the operating part, and the performance cannot be fully exerted.

About Places to Place


Please ensure that there is more than 15cm space behind the massage chair and more than 40cm space in front of the massage chair to avoid touching the wall or other objects when the chair back is tilted.



Attention

- Do not use it in the bathroom or sauna and other places with heavy humidity.
- Do not press the power cord against the base of the massage chair when placing it.
- Do not place the machine in the direct sunlight or near the heating appliances and other high temperature places.
- Do not let open flame objects such as cigarettes or ashtrays near the massage chair:temperature places.
- Do not use it in the place directly facing the air conditioning wind.

- To prevent damage to the floor and reduce running noise, please lay carpets under the massage chair. (Do not lay electric blanket.)
- As the massage chair is heavy, please be careful not to damage the floor when placing it on the ground
- Please place the massage chair 1m away from the TV, radio and other audio-visual electronic equipment. Otherwise, noise may occur
- If a person with a weight of 100kg sits on a massage chair, the weight is about 3.9kg per 1cm

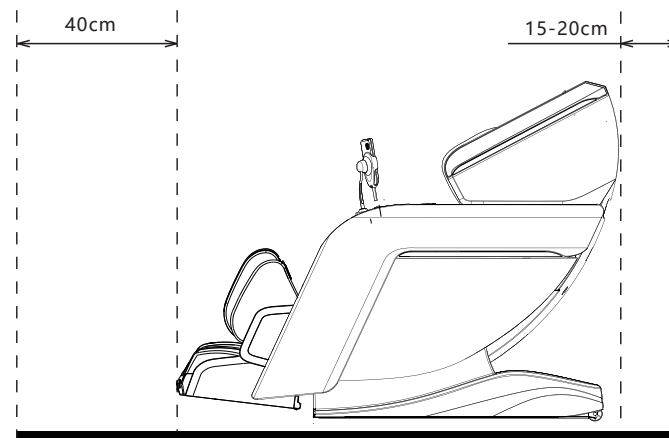


Confirm the items in the packing box

Please confirm that the main body, power cord, instructions and other accessories are installed in the packaging box.

Before Use

Confirm fabric

Open the packing box and confirm whether the cloth of the massage chair is damaged. In addition, please check other parts to ensure that the fabric is not damaged.

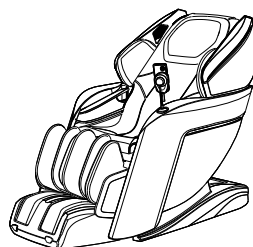


Warning

- Before massage, be sure to confirm whether the fabric at the massage part is cracked or damaged. Please check other parts to ensure that the fabric is not damaged. If you find any damage on the fabric, regardless of the size of the damage, please unplug the power plug and contact the dealer who purchased the product. Please never remove the fabric of the massage area or use it when the fabric is damaged, otherwise there will be serious potential safety hazards.

Method of connecting power

- ① Insert the connection port of the power cord into the access port.
- ② Insert the power plug into the socket.
- ③ Turn on the main power switch behind the massage chair.

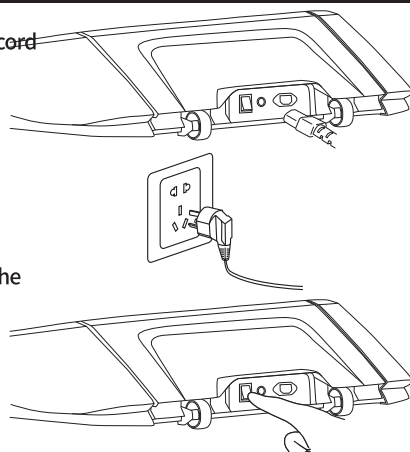

Warning

- Do not use transformers
→ Failure or electric shock may occur.
- Do not use it outside the country where you purchased it.
→ Failure or electric shock may occur.



Warning

- Before inserting the power plug into the socket, make sure that the main power switch is in "O" state.
- Please insert the power plug fully.
- Before use, please follow this instruction manual, turn on all switches in turn, and confirm the normal operation of the machine.

Trouble shooting

Operating sound

The following operating sounds will be made during operation.

With the sound of hammering	Glug glug
Rubbing sound between massage wheel and fabric during kneading operation	Whoosh Squeak
The sound of massage wheel lifting, kneading and finger pressing	Buzz
The sound of air pressure massage	Clunk clunk
The sound of the charging pump	Bang Bang
Airbag exhaust sound	Whoosh
The sound of the chair back tilting up and down	Buzz

※ The above sounds will not affect the performance and down

Safety check

When checking the product under long-term use

Whether the following situations occur

- Sometimes after pressing the switch, it still cannot operate.
- The power cord is cracked or damaged.
- The power cord and power plug are abnormally hot.
- There is a smell of burning or abnormal sound and vibration during operation.
- Sometimes the timer cannot be turned off automatically.
- Other abnormalities or faults occur.

Tips

Love check



If the situation described on the left side occurs, please stop using, unplug the power plug from the socket and send it to the dealer who purchased the product (recorded on the back cover) for inspection to avoid faults or accidents.

For your safe use

The fabric of the massage part is consumables with a service life of 3 years. (Use conditions: 60 minutes per day.) For safety, it is recommended to contact the dealer who purchased this product before the service life is exceeded after purchase to check. If the service conditions are inconsistent with the above, the service life will also change. For safety, it is recommended to conduct a comprehensive inspection of the massage chair after 5 years of purchase. Please consult with the dealer when purchasing this product for the cost of inspection and maintenance.

Q Can the massage chair be placed on the electric carpet?

A It may cause fire, so please do not put it on the electric blanket.
As the massage chair is heavy, it may damage the heating element of the electric blanket, which may cause a fire.

Q Can it be placed on wood floor and floor heating floor?

A Can be placed.
However, the floor may be damaged, so it is recommended to place the massage chair on the carpet. (The same goes for tatami)

Q How to enhance the massage intensity?

A Please try the following operation.
● Fully sit down and place the weight on the back of the chair.
● Adjust the massage force by pressing the force button of the controller.
Through the above methods, the massage intensity can be increased.

Q How to reduce the massage intensity?

A Please try the following operation.
● Adjust the massage force by pressing the force button of the controller.

Q What should I do if I spill water accidentally?

A Electric shock may be caused. Please stop using immediately and do the following.
● Turn off the power switch behind the massage chair
● Pull the power plug out of the socket

After the above operations, please contact the sales shop or dealer (recorded on the back cover) when purchasing this product

Q Will the back rest recline automatically when massage finishes?

A It will not reset automatically when massage finishes if power is turned off
However, when the power is turned off or the timer preventing overuse is turned on, the back rest will reset automatically

Q Can people whose height is less than 150cm or more than 180cm also use it?


A Can be used.

Q Can be use too heavy?

A It is recommended for people weighing less than 100kg

Cleaning, Maintenance and Storage Methods


Cleaning and maintenance methods

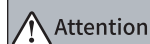

● Be sure to unplug the power plug from the socket before cleaning and maintenance.
In addition, do not plug or unplug the power plug with wet hands

[Cleaning and maintenance of sewing products]

When stained with dirt, immerse the soft cleaning cloth in warm water with a small amount of neutral detergent, screw it completely and then wipe it.

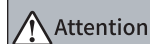

● Discoloration may occur after long contact with plastic and indoor wallpaper.

[Cleaning and maintenance of plastic parts]

Immerse the soft cleaning cloth in warm water containing a small amount of neutral detergent, and completely screw it dry before wiping. To avoid the residue of lotion, please use a dry cloth to wipe it clean again.

● Do not use alcohol, other solvents and polishing powder for cleaning and maintenance, otherwise it may cause scratches, discoloration and cracking

Applicable goods	Not applicable goods

[Cleaning and maintenance of wooden panel, wooden pole and remote control]

Please use a dry cloth to wipe away the dirt on the wooden panel, wooden pole and remote control.

● Do not use a wet towel to wipe.


→ Otherwise, failure may occur.

Storage method

● After decontamination and dust removal, place it in a dry place.

● When not in use, please cover the dust cover to avoid dust.

● To prevent pets from biting the power cord, do not let pets near the machine during storage

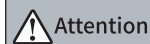

● Do not place the machine in direct sunlight or near heating appliances and other high temperature places.

After-sale service

1. Warranty card


Please be sure to save the product warranty card provided by the sales shop or dealer at the time of purchase

Please read the warranty carefully and keep it properly. (Show it when repairing.) The warranty period is recorded in the warranty card. However, this product is only used in the country where you purchased it, so when it is used abroad, the warranty card is still invalid even within the warranty period. In addition, when it is used for other purposes than general household use, even within the warranty period, the relevant repair costs shall also be paid.

2. When the repair is entrusted

Before commissioning maintenance, if there is still any abnormality after inspection, please contact the sales shop where you purchased this product to commission the following matters, or pay attention to the wechat official account and contact information of the production company on the product warranty card

- 1) Name, address, telephone number
- 2) Product name and model
- 3) For the production number, please refer to the company label affixed on the chair foot cover.
- 4) Purchase date
- 5) Fault status (please inform in detail as much as possible.)

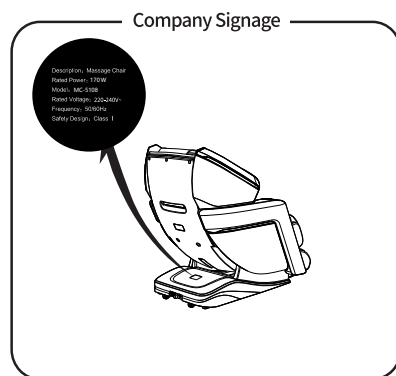

● Within the warranty period

Please show the warranty card during maintenance. Please read the warranty card carefully because it may be a paid repair even within the warranty period.

● Beyond the warranty period expires

Please contact the sales shop when purchasing this product. If the product performance can be maintained through maintenance, paid maintenance can be carried out according to customer needs.

● When it is impossible to entrust the sales shop

Please contact your dealer.

● If you have any questions about after-sales service

Please consult the sales shop or dealer who purchased this product.

● When the product is used in a country other than the country where it was purchased, it may be refused to repair

Specifications

Model	MC-5108
Use power	220-240V~
Rated power consumption	170W
Rated frequency	50-60Hz
Rated time	10~30min
Kneading speed	Automatic program control
Beating speed	Automatic program control
Up and down movement speed of massage wheel	5.0cm/s
Massage range	1010mm
Massage wheel output	10cm
Barometric pressure	About 37kpa
Heater	Below 45 °C
Optional program	20 kinds
Manual program	Multiple coordination
Tilt angle	128 °~156 ° relative to the horizontal plane No grade (electric tilt)
Size of massage chair	About 143*75*117cm(LxWxH)
Weight	About 92.8 kgs
Enclosure	Power cord, manual and other accessories
packing	1 package/ 1set

Children's Safety Notification

⚠ Attention

The massage effect of the massage chair can only be reflected in the normal use of adults. It is recommended that children not use the device. Attention must be paid to the following aspects, otherwise accidents or personal injuries may be caused.

1 Do not let children play on the equipment


Limbs may get stuck by equipment such as



The device may tip over.



When the back of the massage chair is tilted, Children may get stuck in the backrest and floor.

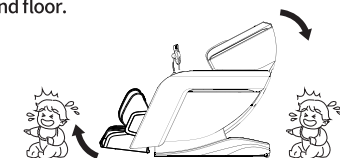

2 To prevent children from using the device alone, slide the power switch to the "O" (off) position when the machine is not in use.

